

VERSANT Sample Preparation 1.0 Reagents Kit



A flexible and reliable solution for extracting viral RNA for COVID-19 testing.

siemens-healthineers.com/versant-sample-prep1



The VERSANT® Sample Preparation 1.0 Reagents Kit* help improve the quality of your nucleic acids extracted from clinical specimens while increasing your lab's efficiency to better manage increasing test volumes. This may be especially critical during the seasonal peak or during a global outbreak such as the COVID-19 pandemic. High-quality nucleic acids result in improved performance of molecular applications that enable you to provide physicians with more reliable test results for patient management.

Extraction of high-quality nucleic acids

- Extraction from routinely used clinical specimens including respiratory samples†
- Nasal swabs, nasopharyngeal swabs, oropharyngeal swabs, nasopharyngeal wash/aspirate or nasal aspirate
- Bronchoalveolar lavage (BAL)



*CE-IVD marked for diagnostic use. Product availability will vary by country.

†See Instructions for Use as to any specific instructions or limitations.


The VERSANT Sample Preparation 1.0 Reagents Kit provides:

In vitro diagnostic (IVD) reagents for clinical use


- Ready-to-use reagents require no preparation
- Shipped and stored at ambient temperature or at 2–8°C (Proteinase K)

Proven technology

- Proprietary magnetic particles coated with a nano layer of silica
 - Results in the efficient isolation of high-quality nucleic acids
- Homogeneous size and shape of magnetic particles
 - Improves reproducibility, recovery, and the quality of nucleic acid
 - Standardizes and enhances performance of molecular applications such as real-time PCR.

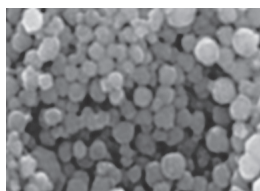


Siemens Healthineers magnetic nano beads

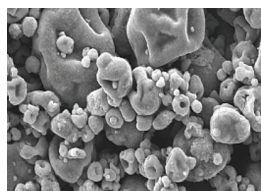

Other commercially available beads

Flexible workflow options

- Automated workflow on the VERSANT kPCR Molecular System with up to 192 tests per workday
- Manual protocol for lower throughput needs

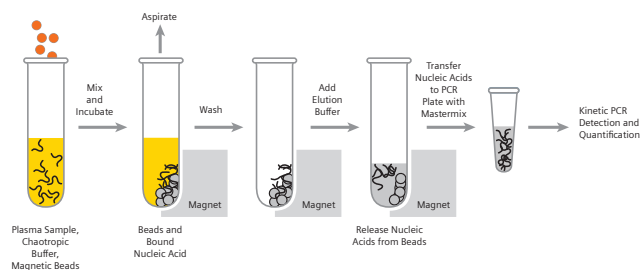

Ordering Information

SMN	Description
10286026	VERSANT Sample Preparation 1.0 Reagents, Box 1 (96 extractions)
10286027	VERSANT Sample Preparation 1.0 Reagents, Box 2 (96 extractions)
10282928	VERSANT kPCR Molecular System SP
10720881	VERSANT kPCR Molecular System AD

Learn more:

siemens-healthineers.com/versant-sample-prep1

VERSANT and all associated marks are trademarks of Siemens Healthcare Diagnostics Inc., or its affiliates. All other trademarks and brands are the property of their respective owners.

Product availability may vary from country to country and is subject to varying regulatory requirements. Please contact your local representative for availability.

Siemens Healthineers Headquarters

Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen, Germany
Phone: +49 9131 84-0
siemens-healthineers.com

Published by

Siemens Healthcare Diagnostics Inc.
Molecular Diagnostics
725 Potter Street
Berkeley, CA 94710-2722
USA
Phone: +1 510-982-4000